



Important! This quick reference guide does not replace the user manual which can be downloaded at www.humancaregroup.com.
Read all instructions before use!

Components

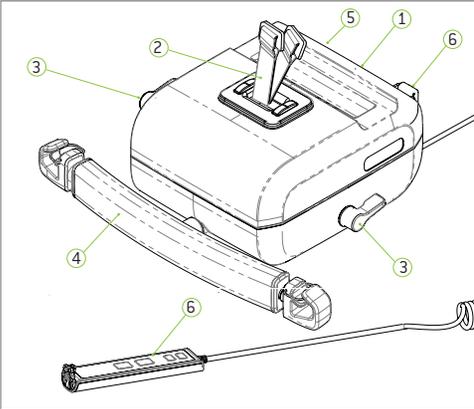

- | | |
|--------------------|--------------------------|
| 1 Secondary Panel | 4 Telescoping Hanger Bar |
| 2 Lift Belts | 5 Emergency Stop |
| 3 Direction Levers | 6 Hand Control |

Secondary Panel Functions

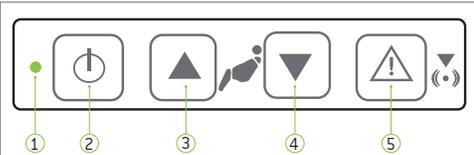


- | | |
|--------------------------|----------------------|
| 1 Status Indicator Light | 4 Lower Lift |
| 2 Power On/Off | 5 Emergency Lowering |
| 3 Raise Lift | |

Status Indicator Lights

- | | |
|--|---|
| Flashing Green
On/Ready for Use
Full Battery | Constant Red
Safety Catch Engaged
Troubleshooting! |
| Flashing Yellow
On/Ready for Use
Charge Battery Warning | Constant Yellow
Lift Charging
Functions Blocked |
| Flashing Red
On/Ready for Use
Low Battery | Flashing Green/Yellow
Service Warning
Perform Inspection |
| Constant Green
Lift in Use
Lift Loaded | Flashing Red/Yellow
Service Passed
Battery Low |

Important! This overhead lift has been designed and tested for indoor use with one patient at a time. This product must not be immersed in water.
Read all instructions for full product information!

Hand Control Functions

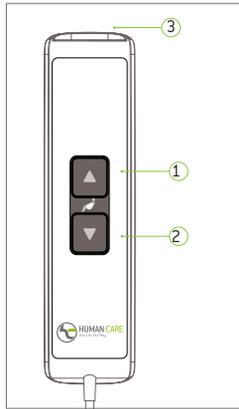


- 1 Up Button
- 2 Down Button
- 3 USB and Charge Port



Important!

The covers for the USB input and the charging input [3] must be closed for the hand control's water protection (IPX5) to be valid.

Using the lift

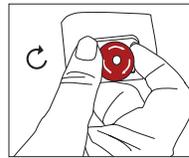

Emergency Stop Activation
Press down on red button to activate.


Emergency De-Activation
Turn red button clockwise.



Important! Emergency lowering must only be used in emergency situations. There will be a 3 second delay.




To activate, press handset "Down" button [2] and secondary panel emergency lowering button [5] at the same time.

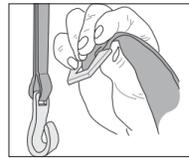

Suspending the lift

Place the lift on a suitable surface and feed out belt by gently pulling out belt while at the same time pressing the "Down" button on either the hand control or secondary panel. Hook lift belt onto distance strap so that it is held firmly in position (easiest to hook from above at an angle as shown). **Ensure the lift belt is firmly in place before use.**

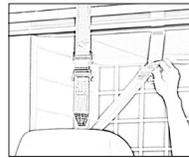


Room to room transfers

Activate desired belt by turning directional lever. If lever faces towards the leg hook, the belt nearest to the leg hook is active. If it faces towards the hanger bar, the belt closest to the bar is active.



Important! Do not turn the lever unless both belts are connected to their respective rail rollers. The lever must always point to the connected belt.

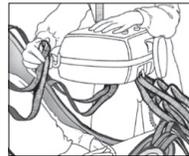

Slings

Attach head sling loops to carry bar first, then leg straps to leg hooks. Ensure the straps are secured beneath the safety pin.



Proper lifting techniques

Do not use the telescoping bar to move the patient and the lift. Do not hold onto the hooks while performing a lift. Once the patient is comfortably secured in sling, raise lift, hold sling and position/move patient.

General Care



Charge lift every night.



Wipe down with slightly dampened drying cloth.



Disinfect with DuPont™ Virkon®.



Never use strong detergents on lift.



Do not immerse in water.



Activate emergency stop when not in use.



Important! Do not use strap to transport!
Lift must always be stored or moved in horizontal position. Rest on cart to move.

Inspection

Annual Inspection - The lift must be inspected at least once per year by a qualified technician certified by Human Care.
Monthly Inspection - Check lift belts for damage or wear. Feed entire length to ensure proper functioning.
Daily Inspection - Always inspect the lift for visible signs of damage before use. Ensure the lift is correctly mounted on the distance strap or rail roller. Always check slings and belts for signs of wear before use. Do not use if wear evident.

Troubleshooting



If the indicator light is not on, check the following:

- Mains switch is on
- Batteries do not require charging
- Emergency stop button is not pushed in (deactivate by turning clockwise)



If the indicator light is flashing, check the following:

- Check the hand control button functions up / down when depressed (replace if necessary)
- When pressing up, the lift is not already in its maximum upper position
- When pressing down, the lift is not already in its maximum lower position



If the indicator light is constant red:

- The lift motor is overheated - allow to cool
- The lift is overloaded - lower using emergency lowering button



If the indicator light is constant yellow:

- Charging in progress



If the lift belt does not feed out properly:

- Pull the belt upwards at the same time as you depress the Down button on the hand control / lift panel




If the indicator changes to flashing red:

- The batteries need charging, patient may be lifted a few more times

Charging the Lift



Flashing Red - Charge lift immediately!

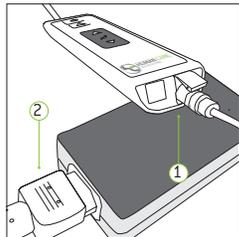

- 1 Connect charger plug to handset.
 - 2 Connect plug to charger and insert plug into wall to charge.
- The charger must only be connected to an earthed socket.
Note: If the emergency stop is activated, the lift will not charge!



Important! The lift can only be used with a Human Care charger and must be disconnected to use the lift.

Models: 55070
Human Care HC AB. All Rights Reserved.

