



The Winco Recliner Collection  
Verō Care Cliner



\*6Y58 with options shown

Designed by Winco, built by you. That's just smart.

Realities of patient care continuously create opportunities for smarter solutions. The **Verō** family of medical recliners has been created to provide total flexibility in matching your application with exactly the features you require and the comfort your patients demand.

The **Verō** recliner's modular design helps the clinician choose a configuration that is perfectly suited to their specific needs:


- Promotes patient independence with low effort recline and infinite position backrest.
- Supports clinician comfort with bilateral, foot-operated controls for most functions.
- Enhances facility efficiency with attractive, durable construction that is robust and easy to clean.

Winco designed it, you build it, that's smart flexibility for smart patient care.



The Winco Recliner Collection  
**Verō Care Cliner**

**ANATOMY OF THE VERō CARE CLINER**

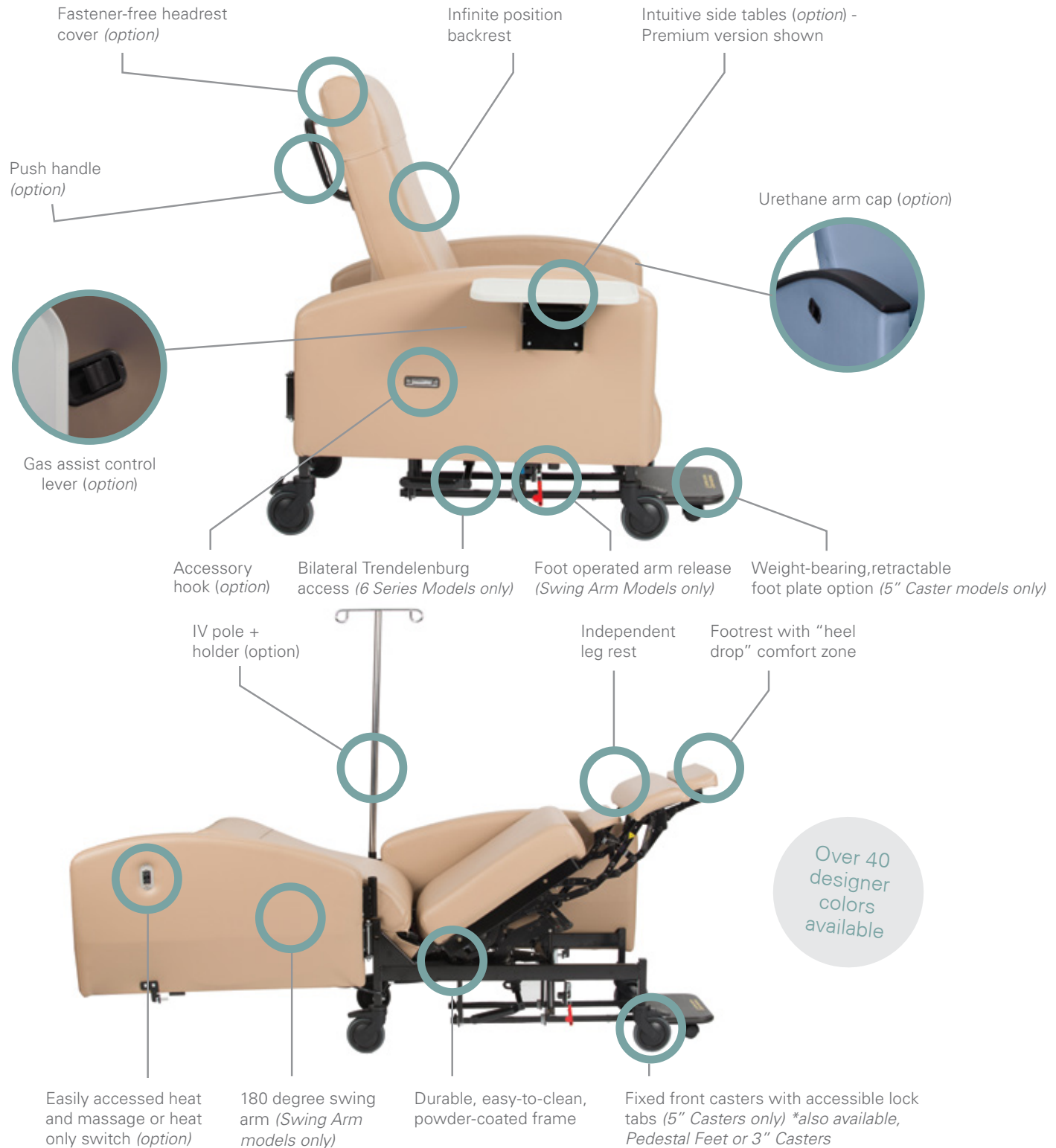

## TRENDELENBURG

Key to clinical utility is the ability to quickly select Trendelenburg positioning where needed. Trendelenburg on **Verō** is controlled with bilateral access foot levers, eliminating the need for the clinician to bend down to reach controls.

Exclusive to **Verō**, Winco introduces "Secure-Glide™" Trendelenburg control. Secure-Glide™ provides gentle patient positioning while preventing clinician strain. The controlled descent helps eliminate any need, or natural reflex, to strain to control the rapid patient positioning in emergency situations.



Secure-Glide™ for stress-free positioning



Foot lever for exclusive Secure-Glide™ Trendelenburg descent

## RECLINE CONTROL OPTIONS

All **Verō** models feature a push back action in order to raise the footrest, and both offer convenient and comfortable infinite positioning of the back rest. **Verō** offers two methods of control of that positioning.



Gas assist shown

**Gas Assist Models** - Our premium gas-back model allows positive locking position of the backrest angle with hydraulic cylinder assisted operation. Locking action is controlled by easily accessed bilateral control levers. Push to extend the footrest and initiate recline, then pull either lever until the desired back angle is achieved. Upon release of the lever, the seat back locks securely in place.

**Push Back Models** - Flowing from the same motion used to extend the footrest, the push back model allows seamless positioning of the desired back rest angle. The gravity operated, gentle resistance designed into the chair then maintains selected positioning. Push back models feature a lower price point and a slightly less encumbered aesthetic.

## ARM OPTIONS

Facility utility is enhanced through the availability of two arm configurations. Standard fixed arms are cost effective, sturdy and durable. Optional swing arms allow 180 degree swing-away movement of the arm for convenient cleaning access and barrier-free patient transfer.

Upholstered arm shown



### Swing Arm Option

Swing arms are ideal when easy access cleaning and transfer of patients is an important consideration.

Each side of the **Verō** with swing arms is easily opened with the convenient foot controls.

Even if an optional table is left in the up position, there is no interference with the full swing motion nor hindrance with patient in the full recline position for true barrier-free transfer.

### Fixed Arm Option

Fixed arms are suitable for most general purpose operations when transfer of patient is not required or where working space for swing arm motion is severely limited.

### Urethane Cap Option

The Urethane cap option provides unique styling, firm support, and extended arm cap life.



Urethane Arm Cap (option) >>

## MATCH YOUR PURPOSE

The **Verō** modular design allows choices in how the chair is anchored or moved. Whether stationary or for patient transport, **Verō** has a perfect solution. **Verō** moves with the configuration you need.



**Pedestal Feet** - Looking for a stable room recliner for guests, or don't want your residents moving chairs around? Pedestal feet are the solution. When mobility is not required or desired, pedestal feet provide a less clinical look and a lower seat height. Pedestal feet create a seat height of 17.5" from the floor, meeting ADA Compliance.\*

**Casters** are available for ease of movement of the chair for cleaning, maintenance and staging for patient access. All casters come with lockable, foot lever release controls.

*It is recommended that after moving the chair always lock the casters prior to patient ingress/egress.*



Pedestal



3" Casters



5" Casters

**3" Casters** - Want a low seated height but need mobility for cleaning and room configuration flexibility? Then 3" casters are probably best. 3" casters have a seat height of 20.5" from the floor meeting ADA Compliance\*.

**5" Casters** are first choice when moving the chair about a facility. 5" casters have a seat height of 22.75" from the floor. Want a little more height to bring the patient closer to staff? 5" casters are your answer.

\* Reference ADAAG 4.19.2 Heights and Clearances.

## FOOT PLATE

With 5" casters, patients will appreciate the added comfort of the retractable foot plate. Paired with the pushbar option, the **Verō** with footplate is suitable for within-facility patient transport.

**Verō's** footplate is designed to be weight bearing and won't tip when weight is placed on the footplate. The footplate locks as weight is applied to prevent unwanted or unexpected movement.

*Patients should be assisted at all times by a staff member when using the footplate.*



*Footplate is easily positioned by foot.*

## PUSH HANDLE OPTION

The push handle option is recommended with the footplate models if the chair is to be used for in-facility patient transport.



**CLINICAL UTILITY AND FACILITY FLEXIBILITY ARE MAXIMIZED BY VERŌ SIZE/CAPACITY OPTIONS**

**Verō's** innovative design and enhanced utility is available in two sizes: standard and X-large. No matter which size is specified, **Verō** offers comfort, beauty and capability.

- Standard seat width is 22.5" between arms with a 350 pound weight capacity. This chair makes petite to average size patients feel right at home, comforted and individually cradled.
- X-large seat width is 26.5" between arms with a 500 pound weight capacity. This increased size makes it possible to care for higher BMI patients while providing a comfortable and personalized experience.



*X-Large 500 lb Capacity*

*Standard 350 lb Capacity*

**Table Options** There are three table options; standard, pivot and premium. All tables on the **Verō** are available in left, right or both configurations and can be ordered in a variety of combinations to assist patient comfort and clinician preference. All feature a stain resistant, non-porous material and all but the pivot table have a lip edge design to contain small spills.

- **Standard table** This value priced lift-to-raise-and-lower table will be familiar to current Winco Elite series owners.
- **Pivot table** swivels so patient can have the table in front of them.
- **Premium table** features an intuitive single hand up/down feature with audible positive lock and release feedback.



*Standard Table*



*Pivot Table*



*Premium Table*

Standard tables add 2.75" (2.5" for XL) to each side of the chair when lowered, Premium tables add 1.25" (1" for XL) to each side of the chair when lowered. This added width should be considered when ordering table options appropriate for your facility.

**COMPLEMENT DESIGN AND UTILITY BY ADDING THE APPROPRIATE ACCESSORIES**

- **Heat option** provides patient-selected, soothing warmth for maximum comfort.
- **Heat and massage** adds a selection of gentle vibrating massage modes to the heat option, enhancing patient comfort level.
- **Headrest cover** is standard on Trendelenburg models, and is available for all other models as an option. Covers extend the life of the headrest, and are fastener-free, making them easy to lift on and off and keep the chair clean.
- **Accessory hook** is available for either left, right or both sides. The accessory hook is convenient for uses such as a Foley bag. The marine grade stainless steel remains flush until needed to prevent snagging other objects.
- **IV Pole and holder option** for either left, right or both sides. The holder moves with the swing arm, without interfering with chair operation.



*Conveniently placed controls for patient independence*

*Heat & Massage: 3 levels of heat, 5 levels of massage*

*Heat Only: Low/Med/High*



*Accessory Hook (shown here with Foley bag)*



*IV holder mount allows IV pole to remain relatively stationary when the swing arm (option) is in operation*

**MODELS & SPECIFICATIONS** \* Multiple options are available. Please contact us to ask about specifications for your exact configurations.

	PEDESTAL FEET		3" CASTERS		5" CASTERS		FOOTPLATE (5" Casters Only)	
	STANDARD	X-LARGE	STANDARD	X-LARGE	STANDARD	X-LARGE	STANDARD	X-LARGE
<b>MODELS - Push Back</b>	5Y04, 5Y44, 6Y04, 6Y44	5Y14, 5Y54, 6Y14, 6Y54	Y03, 5Y43, 6Y03, 6Y43	5Y13, 5Y53, 6Y13, 6Y53	5Y00, 5Y40, 6Y00, 6Y40	5Y10, 5Y50, 6Y10, 6Y50	5Y08, 5Y48, 6Y08, 6Y48	5Y18, 5Y58, 6Y18, 6Y58
<b>MODELS - Gas Back</b>	5X04, 5X44, 6X04, 6X44	5X14, 5X54, 6X14, 6X54	5X03, 5X43, 6X03, 6X43	5X13, 5X53, 6X13, 6X53	5X00, 5X40, 6X00, 6X40	5X10, 5X50, 6X10, 6X50	5X08, 5X48, 6X08, 6X48	5X18, 5X58, 6X18, 6X58
<b>Seat Height to Floor</b>	19.25" (48.9 cm)	19.25" (48.9 cm)	20.5" (52.1 cm)	20.5" (52.1 cm)	22.75" (57.8 cm)	22.75" (57.8 cm)	22.75" (57.8 cm)	22.75" (57.8 cm)
<b>Seat Height to Footplate</b>	-	-	-	-	-	-	19.5" (49.5 cm)	19.5" (49.5 cm)
<b>Seat Width (Between Arms)</b>	22.5" (57.2 cm)	26.5" (67.3 cm)	22.5" (57.2 cm)	26.5" (67.3 cm)	22.5" (57.2 cm)	26.5" (67.3 cm)	22.5" (57.2 cm)	26.5" (67.3 cm)
<b>Seat Depth</b>	22" (55.9 cm)	22" (55.9 cm)	22" (55.9 cm)	22" (55.9 cm)	22" (55.9 cm)	22" (55.9 cm)	22" (55.9 cm)	22" (55.9 cm)
<b>Overall Height</b>	44.75" (113.7 cm)	44.75" (113.7 cm)	46" (116.8 cm)	46" (116.8 cm)	48.25" (122.6 cm)	48.25" (122.6 cm)	48.25" (122.6 cm)	48.25" (122.6 cm)
<b>Overall Width (w/out table &amp; w/ fixed arm or w/ swing arm closed,)*</b>	29.5"	34"	29.5"	34"	29.5"	34"	29.5"	34"
<b>Width: Swing arm open 90°</b>	+35.75"	+35.75"	+35.75"	+35.75"	+35.75"	+35.75"	+35.75"	+35.75"
<b>Width: Swing arms open 180°</b>	+28.5"	+32.5"	+28.5"	+32.5"	+28.5"	+32.5"	+28.5"	+32.5"
<b>Back Height (Above Seat)</b>	29.25" (74.3 cm)	29.25" (74.3 cm)	29.25" (74.3 cm)	29.25" (74.3 cm)	29.25" (74.3 cm)	29.25" (74.3 cm)	29.25" (74.3 cm)	29.25" (74.3 cm)
<b>Overall Length (Upright)</b>	39.5" (100.3 cm)	39.5" (100.3 cm)	39.5" (100.3 cm)	39.5" (100.3 cm)	39.5" (100.3 cm)	39.5" (100.3 cm)	39.5" (100.3 cm)	39.5" (100.3 cm)
<b>Overall Length (Reclined)</b>	69" (175.3 cm)	69" (175.3 cm)	69" (175.3 cm)	69" (175.3 cm)	69" (175.3 cm)	69" (175.3 cm)	69" (175.3 cm)	69" (175.3 cm)
<b>Arm Height</b>	25.5" (67.3 cm)	25.5" (67.3 cm)	26.75" (67.9 cm)	26.75" (67.9 cm)	29" (73.7 cm)	29" (73.7 cm)	29" (73.7 cm)	29" (73.7 cm)
<b>Weight (no options)*</b>	104 lb (47.2 kg)	122 lb (55.3 kg)	104 lb (47.2 kg)	122 lb (55.3 kg)	104 lb (47.2 kg)	122 lb (55.3 kg)	114 lb (51.7 kg)	132 lb (59.9 kg)
<b>Weight Capacity</b>	350 lb (158.8 kg)	500 lb (226.8 kg)	350 lb (158.8 kg)	500 lb (226.8 kg)	350 lb (158.8 kg)	500 lb (226.8 kg)	350 lb (158.8 kg)	500 lb (226.8 kg)

**FEATURES**

NON-TRENDELEBURG (5 SERIES)	ARM TYPE	PEDESTAL FEET	3" CASTERS	5" CASTERS	FOOTPLATE 5" CASTERS
<b>STANDARD (350 LB CAPACITY)</b>	Fixed Arm	5Y04	5Y03	5Y00	5Y08
	Swing Arm	5Y44	5Y43	5Y40	5Y48
<b>XL (500 LB CAPACITY)</b>	Fixed Arm	5Y14	5Y13	5Y10	5Y18
	Swing Arm	5Y54	5Y53	5Y50	5Y58

TRENDELEBURG (6 SERIES)	ARM TYPE	PEDESTAL FEET	3" CASTERS	5" CASTERS	FOOTPLATE 5" CASTERS
<b>STANDARD (350 LB CAPACITY)</b>	Fixed Arm	6Y04	6Y03	6Y00	6Y08
	Swing Arm	6Y44	6Y43	6Y40	6Y48
<b>XL (500 LB CAPACITY)</b>	Fixed Arm	6Y14	6Y13	6Y10	6Y18
	Swing Arm	6Y54	6Y53	6Y50	6Y58

**OPTIONS**

	TB 133	HEAT	HEAT & MASSAGE	STANDARD TABLES (LIFT & DROP)	PREMIUM TABLES (SQUEEZE & DROP)	PIVOT TABLE	IV POLE	ACCESSORY HOOK	PUSH HANDLE	ARM REST COVERS	HEAD REST COVERS	MEMORY FOAM UPGRADE	URETHANE ARM CAPS
<b>AVAILABLE ON ALL MODELS</b>	All	All	All	All; L, R, Both	All; L, R, Both	All; L, R, Both	All; L, R, Both	All; L, R, Both	Recommended for Footplate models only	All	Standard on 6 Series, Optional on 5 Series	All	All

**WARRANTY** Contact Winco or visit our website for warranty details.

**COLORS** Contact Winco or visit our website for available color selections.

**NOTE** Winco reserves the right to change specifications without notice.

**CERTIFICATIONS** **1)** The Verō Care Cliner meets ANSI/BIFMA X5.4 for reliability testing **2)** Winco foam meets California Technical Bulletin 117, Flammability Test Procedure for Resilient Filling Materials used in Upholstered Furniture **3)** TB 133 Option: Meets flammability standards of California Technical Bulletin 133, Flammability Test Procedure for Seating Furniture for Use in Public Occupancies and ASTM E1537, Standard Test Method for Fire Testing of Upholstered Furniture **4)** Heat and Massage, Heat Only Options: Meets electrical safety standards, UL 1647 and CAN/CSA-C22.2 No.68-92



**Contact Us**

To arrange a demonstration or to speak with one of our associates, please feel free to contact us at:  
**800-237-3377 | WincoMfg.com | 5516 SW 1st Lane, Ocala, FL 34474 | Info@WincoMfg.com**