

# LIFTKAR PTR

POWERED TRACKED STAIRCLIMBER



# LIFTKAR PTR - POWERED TRACKED STAIRCLIMBER

The Liftkar PTR is a tracked stairclimber that helps you safely convey wheelchair users up and down stairs.



Control panel handle unit



Adjustable backrest clamp



Control panel base unit


## STANDARD ACCESSORIES



Lap belt PTR  
Art. No. 075 105  
(fix mounted)



Charger PTR  
Art. No. 075 637

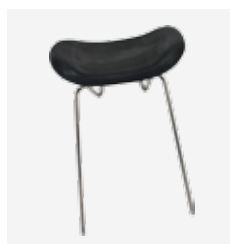

Headrest PTR  
Art. No. 075 104



Backrest cushion PTR  
Art. No. 975 104



Ramp PTR  
Art. No. 975 103

## OPTIONAL ACCESSORIES

## FEATURES:

- Electronic overload protection
- Electronic tilting sensor incl. speed control
- Acoustic and visual driving angle warning
- Microprocessor controlled operation system
- Integrated diagnosis system
- Operation time read out system
- High autonomy: up to 50 floors
- Drive direction LED
- Key switch
- Emergency stop pushbutton




Machine including options and foldable wheelchair-support.

**The LIFTKAR PTR is a mobile unit that has been especially developed for operating on straight steps. Wheelchair users can now move around more freely and helps no longer need to strain their backs on the stairs.**

It's so easy - simply release two levers to dismantle the LIFTKAR PTR into two easily portable parts. The base unit with battery pack fitted weighs approx. 28 kg, while the handle unit weighs approx. 11 kg.

The Liftkar PTR ensures extremely safe handling on the stairs.

Total machine weight:	39 kg
Handle weight:	11 kg
Climbing unit weight:	28 kg
Max load (wheel chair and passenger):	130 kg / 160 kg
3 speeds (per minute):	4.5 m / 5.5 m / 6.5 m
Max. stair incline:	35°
Height:	930 mm
Width:	650 mm
Overall length base unit:	980 mm
Overall length incl. handle unit:	1380 mm
Landing space dimensions (min.):	970 mm x 970 mm

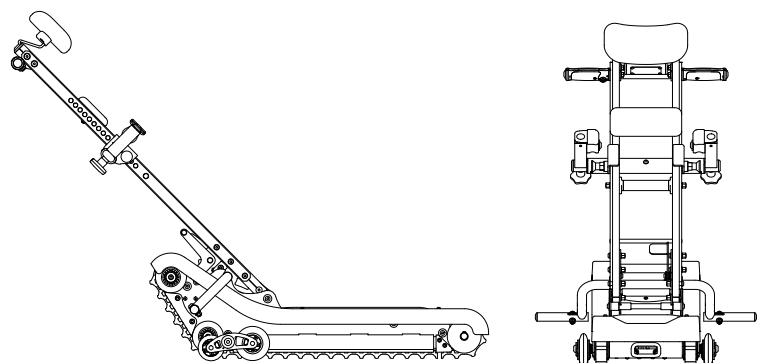



**SANO Transportgeraete GmbH**  
Am Holzpoldlgut 22  
4040 Lichtenberg / Linz  
Austria

Tel. +43 7239 / 510 10  
Fax +43 7239 / 510 10 14  
office@sano.at  
[www.sano.at](http://www.sano.at)

**SANO Deutschland GmbH**  
Bahnhofstrasse 5  
83080 Oberaudorf  
Germany

Tel. +49 8033 / 308 96 0  
Fax +49 8033 / 308 96 17  
info@sano-treppensteiger.de  
[www.sano-deutschland.de](http://www.sano-deutschland.de)

**SANO UK Powered Stairclimbers Ltd.**  
Bristol Court, Betts Avenue  
Martlesham Heath  
Ipswich, Suffolk / IP5 3RY, England

Tel. +44 1473 / 333 889  
Fax +44 1473 / 333 742  
info@liftkardirect.com  
[www.liftkardirect.com](http://www.liftkardirect.com)