


CAMINO

TREATMENT/
BOBATH/
MASSAGE TABLES

mingtenghealth.com
Edition of 2020

Shanghai Sintang Industrial Co.,Ltd
Since 2003
Edition of 2020

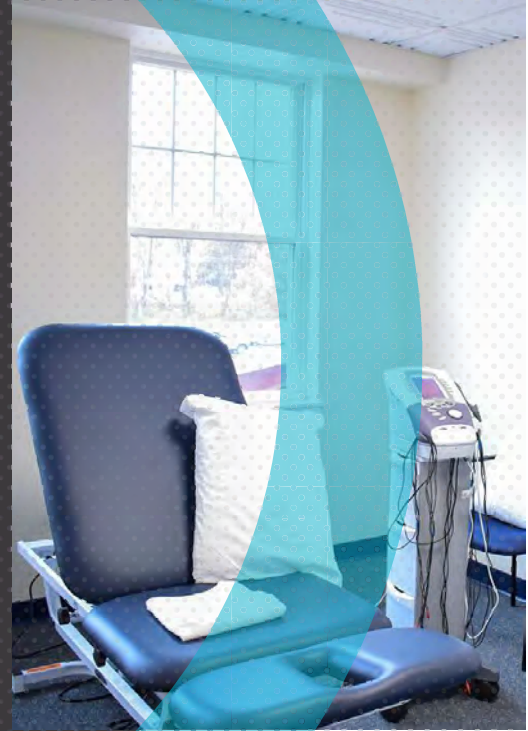

Treatment / Bobath / Massage

Table /
////////





CAMINO

Treatment / Bobath / Massage Tables



Treatment Table Series



Bobath Table Series



Massage Table Series

Today, Camino has 12 different models and more than 10 optional accessories, which can be customized according to the user's habits. It can be used in a broader application scenario. Therefore, in the past 10 years, Camino has become the favorite classic design of masseurs, physiotherapists, chiropractors, physician and surgeons.

Camino Treatment Basic


C	T	B	-									
1	2	3		4	5	6	7	8	9	10	11	

CTB is the first model of Camino series, and it is also a classic design normally used by physiotherapists. It can be used for many purposes such as consulting room examination, massage, B-ultrasound, bone setting, spine treatment, etc. This utility model has a pneumatic adjustable head, which can adjust an angle of $-75^{\circ} \sim 40^{\circ}$, and can be used for leg inspection or as a headrest.



Specification

- Length: 78"/198cm
- Width: 22"/56cm, 25"/64cm, 27"/69cm, 29"/74cm
- Height: 18"/46cm~37"/95cm
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 136lbs/62kg

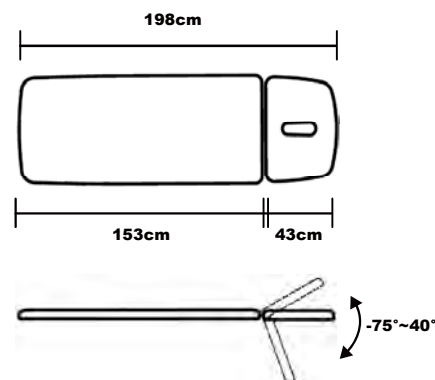

Posture Varieties:




Camino Treatment Infinity

C	T	I	-									
1	2	3	4	5	6	7	8	9	10	11		


The Camino Treatment Infinity multi-functional diagnosis and treatment bed offers more posture adjustments. In addition to the -75° to 40° pneumatic head section adjustment and the 0°~80° tail section angle adjustment, it also has a 40° middle section adjustment. It provides both users and doctors with a comfortable posture regardless of whether the user is in a sitting position, supine position, or prone position.

Specification

- Length: 78"/198cm
- Width: 22"/56cm, 25"/64cm, 27"/69cm, 29"/74cm
- Height: 18"/46cm~37"/95cm
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 174lbs/79kg

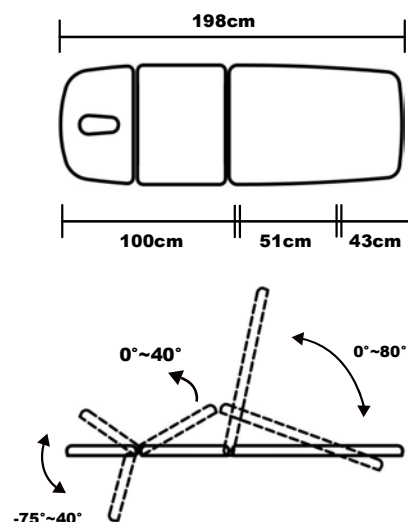

Posture Varieties:

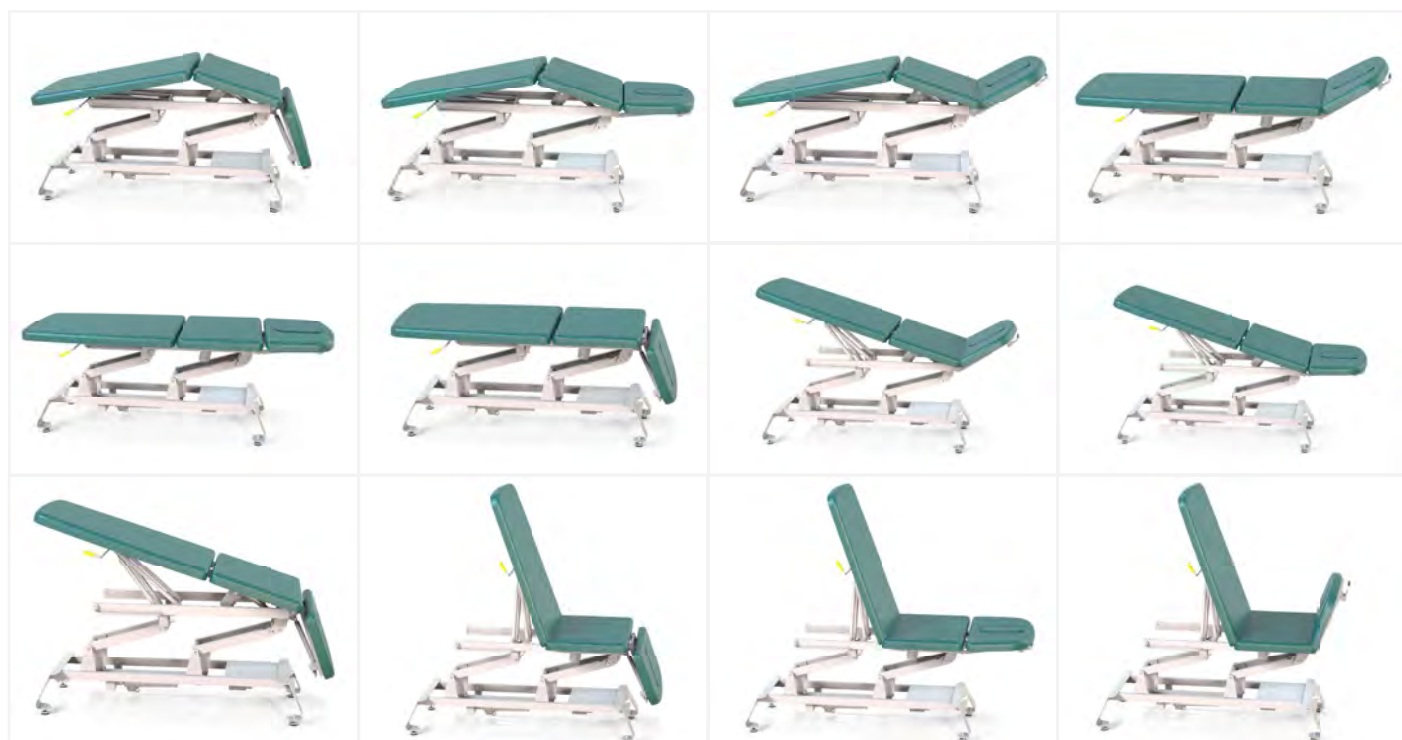

Camino Treatment Scissors


CTS	-											
1	2	3	4	5	6	7	8	9	10	11		

CTS with its Leg First feature can provide better body posture for leg and waist examination and treatment. Two tail section can provide independent pneumatic controlled angle adjustment of $0^{\circ}\sim 80^{\circ}$. Its -75° to 40° pneumatic angle adjustable head piece, can help patients to maintain a good sitting, prone or supine posture.



Specification

- Length: 78"/198cm
- Width: 27"/69cm, 29"/74cm
- Height: 18"/46cm~37"/95cm (Dewert)
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 176lbs/80kg

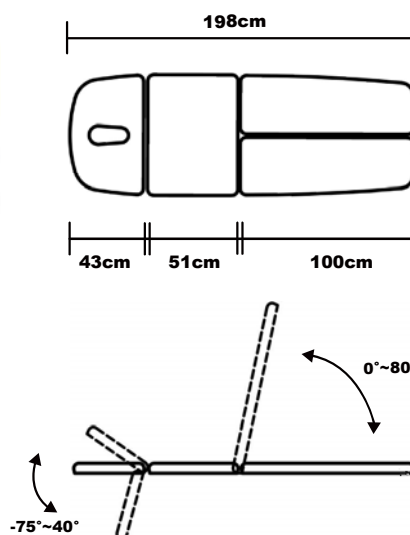

Posture Varieties:



Camino Treatment Russell


C	T	R	-									
1	2	3		4	5	6	7	8	9	10	11	

CTR is a multi-functional diagnosis and treatment bed with high popularity. 50% of the users are from massage and beauty market, and 50% of the users are from medical field. Massage and beauty users tend to use a 3" thick and soft sponge system, while medical users tend to use a 2-2.5" thick hard sponge system. The $-10^{\circ} \sim 80^{\circ}$ tilting backrest angle can meet the needs of supine and prone users.



Specification

- Length: 78"/198cm
- Width: 22"/56cm, 25"/64cm, 27"/69cm, 29"/74cm
- Height: 18"/46cm~37"/95cm
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 147lbs/67kg

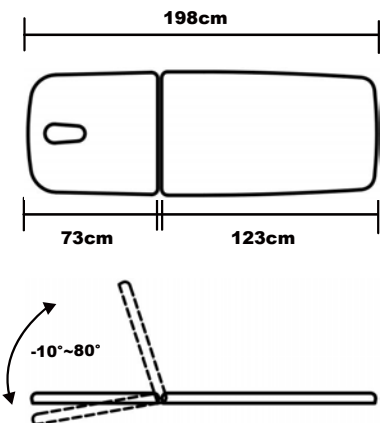

Posture Varieties:



Camino Treatment Danvers


CTD	-											
1	2	3	4	5	6	7	8	9	10	11		

The CTD 3 section treatment / therapy table is a very strong and versatile three-section couch featuring an adjustable angle backrest and foot section support for positioning the patient for various sitting or lying procedures. Both the backrest and foot section are easily adjustable, featuring a smooth anti-entrapment gas operated back & foot rest, and Hemet Medical guarantee this table forms the backbone of many clinics.




Specification

- Length: 78"/198cm
- Width: 22"/56cm, 25"/64cm, 27"/69cm, 29"/74cm
- Height: 18"/46cm~37"/95cm (Dewert)
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 147lbs/67kg

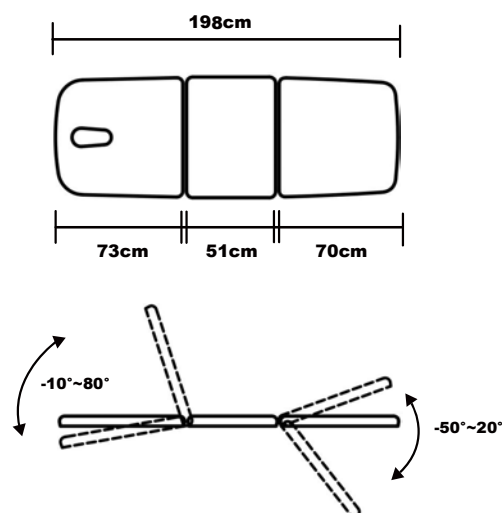

Posture Varieties:

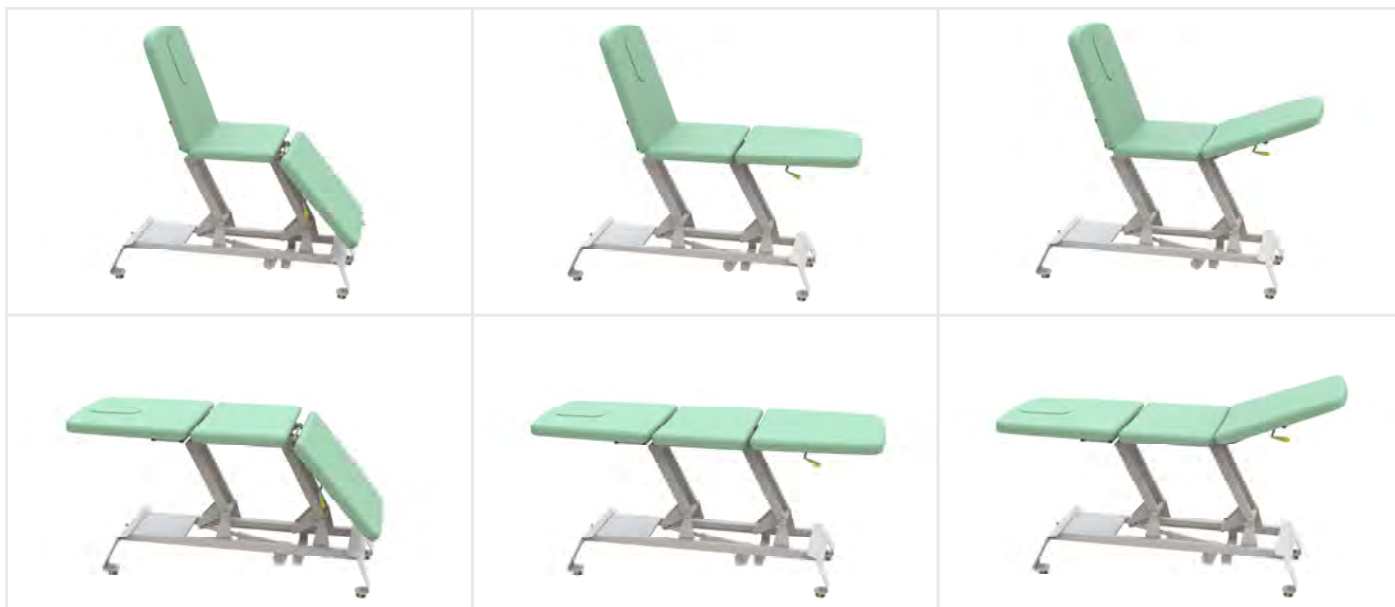

Camino Treatment Avalon



C	T	A	-									
1	2	3	4	5	6	7	8	9	10	11		

This CTA electric treatment table with 3 sections and 2 electric motors allows for many configurations. The central fold, adjustable by electric motor, provides great versatility and many possibilities. This is a very strong and rigid bed designed for various manipulation and mobilization procedures and techniques. This model features a pneumatically adjusted backrest and an electrically adjusted middle section, providing additional comfort for patients.



Specification

- Length: 78"/198cm
- Width: 22"/56cm, 25"/64cm, 27"/69cm, 29"/74cm
- Height: 18"/46cm~37"/95cm (Dewert)
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 167lbs/76kg

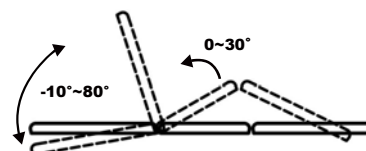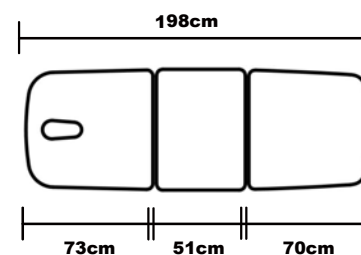

Posture Varieties:



Camino Bobath Flat

C	B	F	-								
1	2	3		4	5	6	7	8	9	10	11

Designed for treatment of adults and children with diagnosed problems with the central nervous system, motor neuron disorders and patients with loss of sensory ability or poor balance. It is used also for Bobath and Vojta therapy, orthopedic treatment or therapy after apoplexy. The generous bed space gives patients and therapists considerable room to complete various rehabilitation exercises.



Specification

- Length: 79"/200cm
- Width: 39.5"/100cm, 47"/120cm
- Height: 18"/46cm~37"/95cm (Dewert)
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 176lbs/80kg

Camino Bobath Tilt


C	B	T	-								
1	2	3		4	5	6	7	8	9	10	11

Specially designed for treating patients with neurological conditions. With a treatment area of 120cm x 200cm and low 46cm minimum height provides excellent access for patients with limited mobility or wheelchair transfers. The electric lifting bed comprises a lifting bed body and an adjustable backrest disposed on the bed body. The backrest section positioned by imported double structured synchronous pneumatic springs, safe and reliable.



Specification

- Length: 79"/200cm
- Width: 39.5"/100cm, 47"/120cm
- Height: 18"/46cm~37"/95cm (Dewert)
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 185lbs/84kg

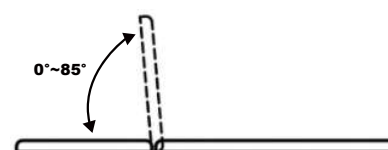

Bobath Round Bolster



Bobath Round Bolster

Bobath bolster adopts hard styrofoam round mandrel, wrapped with 2cm medium hardness recycled sponge and PU vinyl cover on the surface. It can be used for balance and coordination training in hemiplegia, cerebral palsy and other dyskinesia patients

Φ20cm * 80cm (Model: BR2080)

Φ20cm * 160cm (Model: BR2016)

Φ30cm * 80cm (Model: BR3080)

Φ30cm * 160cm (Model: BR3016)

Φ40cm * 80cm (Model: BR4080)

Φ40cm * 160cm (Model: BR4016)

Camino Leg Support Bracket

Leg Support Bracket

(Model Code: CSB01)

Leg support bracket can be installed on the Bobath table to provide patients with better leg and body support. The height-adjustable table mat can be adjusted according to the needs of treatment.



Tips: Tips: Leg Support Bracket can be installed on all sizes of Camino Flat and Camino Tilt.

POSITION CODE 4

Control System

- A - Single Pedal
- B - Double Pedal
- C - Single Pedal+ Side Control
- D - Single Pedal+ 2 Button Hand Control
- E - Single Pedal+ 4 Button Hand Control
- F - Hand Control
- G - Walk-Free Control
- H - Walk-Free Control + Side Control
- I - Walk-Free Control + Hand Control

POSITION CODE 5

Breath Hole Design

- A - No Breath Hole
- B - One Breath Hole + Filler
- C - One Breath Hole + Neck Pillow
- D - Two Breath Hole + Filler
- E - Face Cushion Cut Out Design
- F - Two Breath Hole + Neck Pillow

POSITION CODE 6

Front Arms

- 1 - With Height Adjustable Front Arms
- 2 - Rotatable and Adjustable Front Arms
- 0 - Without Front Arms

POSITION CODE 7

Side Armrest

- 1 - With Revolving Side Armrest
- 2 - With Revolving Guardrail
- 3 - With Rotatable Side Armrest
- 0 - Without Anything

POSITION CODE 8

Stand-by Castor System

- 1 - With Stand-by Castor System
- 0 - Without Stand-by Castor System

POSITION CODE 9

Electrical Safety Lock

- 1 - With Electrical Safety Lock
- 0 - Without Electrical Safety Lock

POSITION CODE 10

Paper Roll Holder

- 1 - With Paper Roll Holder
- 0 - Without Paper Roll Holder

POSITION CODE 11

Sponge System

- 1 - Soft Polyurethane Foam
- 2 - Hard Polyurethane Foam
- 3 - Self-Skin Injection Foam

Model:

C	T	C	-	G	B	1	1	1	1	1	3
---	---	---	---	---	---	---	---	---	---	---	---


Code Position:

1	2	3	-	4	5	6	7	8	9	10	11
---	---	---	---	---	---	---	---	---	---	----	----

Model Codes

Option Codes

Camino's product code consists of three model codes, eight option codes. The seven option codes for Camino table depend on the configuration of various options available. Each position in the code has a specific meaning, and corresponds to one or more characteristics of the table. Our standard model of Camino itself does not contain Walk-free Control, Front Armrests, Side Armrests, Guardrails, Stand-by Castor System, Paper Roll Holder, Electrical Safety Lock, all optional parts need to be purchased separately. All standard models adopt polyurethane sponge mattress, and self-skin sponge system is allowed for limited models and sizes.

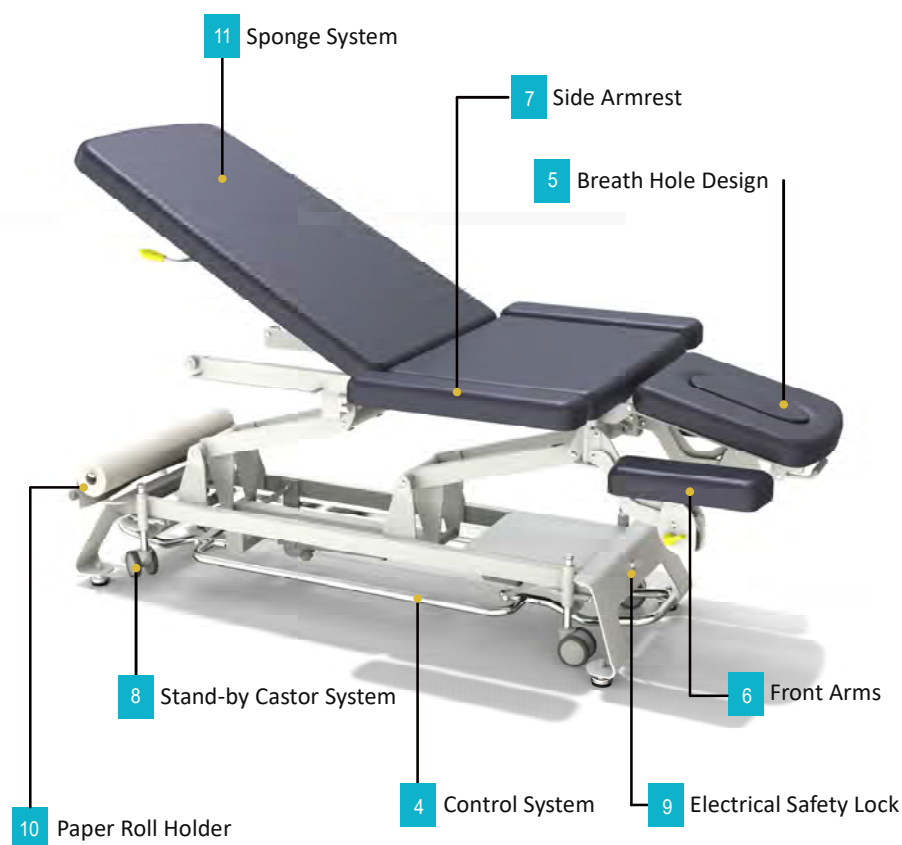

Revolving Guardrail




Face Cushion Cut Out & Filler



Electrical Safety Lock

POSITION 4

Control System

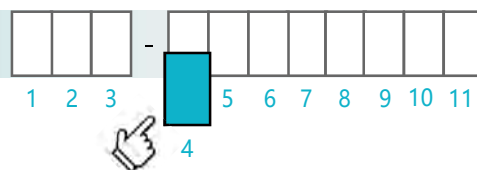

We provide a variety of motor control options, all controllers must be installed in advance by original factory. Please indicate in advance when ordering.



A Single Pedal

Standard controller of CTB, CTC, CTS, CTR, CTD, CBF, CBT, CMF, CMT. It is used to raise or lower the height of the table.



B Double Pedal

Standard controller of CTI & CTA, CMD. It's responsible to raise/lower the height of the table and control the middle section's lifting angle.



C Single Pedal + Side Control

Foot Pedal is for adjusting the height of the table, Side Control is for raising or lowering the middle section of CTI, CTA & CMD.



D Single Pedal + 2 Button Hand Control

It is applicable to CTB, CTC, CTS, CTR, CTD, CBF, CBT models. Foot Pedal and Hand Control are used to adjust the height of the bed. Not available for CTI, CTA & CMD.




E Single Pedal + 4 Button Hand Control

For CTI, CTA & CMD, Foot Pedal is for adjusting the height of the table, hand control is for adjusting the height of the table as well as the mid-lift of the table.



F Hand Control

2 Button Hand Control is for adjusting the height of the table. 4 Button Hand Control is for raising/lowering the height of the table and controlling the mid-lift of the table.

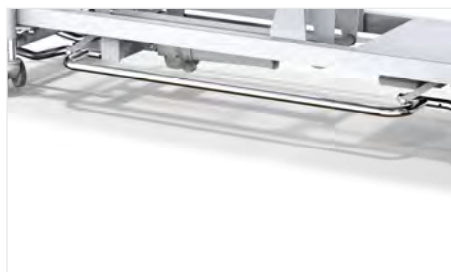


G Walk-Free Control

Walk-Free Control makes height adjustment easily from all the way around the bed. It is applicable to CTB, CTC, CTS, CTR, CTD, CBF, CBT models.



H Walk-Free Control + Side Control

Walk-Free Control is for adjusting the height of the table. Side Control is for raising or lowering the middle section. It is applicable to CTI, CTA & CMD models.

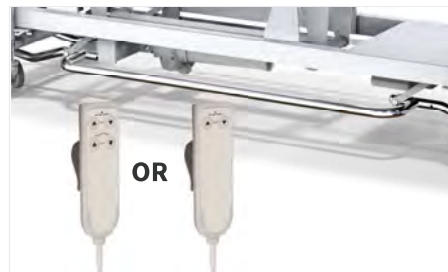

I Walk-Free Control + Hand Control

Walk-Free Control is for adjusting the table's height. 4 Button Hand Control is for raising/lowering the height of the table and the middle section of the table.

POSITION 5

1	2	3	4	5	6	7	8	9	10	11	

Breath Hole Design

The breath hole and face cushion cut-out design were set in advance by original factory. Please indicate in advance when ordering.



A No Breath-hole



B 1 Breath Hole + Filler



C 1 Breath Hole + Neck Pillow

It is applicable to all treatment beds without front arms, or treatment beds with front arms but 29" wide. (Bobath not applicable)



D 2 Breath Hole + Filler

⚠ Suitable for Bobath table only!



E Face Cushion Cut-out + Face Cushion + Filler

⚠ Suitable for Basic, Cabell, Infinity, Scissors



F 2 Breath Hole + Neck Pillow

⚠ Suitable for Bobath table only!

POSITION 6

1	2	3	4	5	6	7	8	9	10	11

Front Armrests

The height adjustable front armrests can adjust the height in the range of 0-20cm through pneumatic pump control, and can be locked at any position. It is used to support the arms of the patient when patient is in prone position. The rotatable & height adjustable front arms are fully adjustable in height and position since the arm frame can rotate and the arm cushion can rotate by itself.



1 With Height Adjustable Front Arms



2 With Rotatable & Height Adjustable Front Arms




0 Without Front Arms

Tips: Tips: Front Armrest is only suitable for Camino Treatment Beds with width of 27" or more.

Rotatable & Height Adjustable Front Arms are not applied on Treatment Infinity, Bobath Flat and Bobath Tilt.

For Treatment Basic, Cabell, Scissors tables, Rotatable & Height Adjustable Front Arms cannot be installed on a table using Revolving Guardrail or Built-in Rotating Side Armrest.

POSITION 7

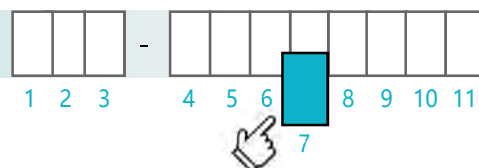


Side Armrest & Guardrails




1 With Revolving Extended Side Armrest (Part code: SA02)

When extended side armrests are rotated to horizontal position, it will provide patient larger support width and effectively improve the comfort and security of users. When you need to get closer to patient for treatment, you can flip armrest down. Side Armrest can be factory installed or purchased separately.

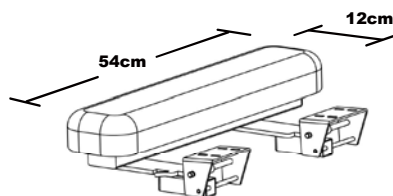


2 With Revolving Guardrail (Part code: RG04)

Warning: Cannot be installed on a table using Rotatable & Height Adjustable Front Arms. Guardrail system will prevent patient from falling or rolling out of the bed, help patient reposition, sit-up, balance, stand or transfer to a table. It can be installed in factory or by customer. The swing down feature allows this guardrail to pivot down and out of the way when getting in and out of bed or removing sheets.

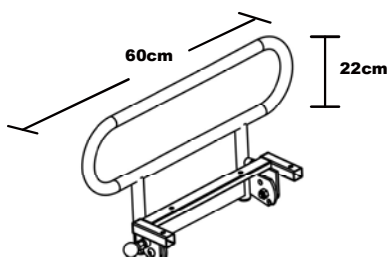


3 With Built-in Rotating Side Armrest

Warning: Cannot be installed on a table using Rotatable & Height Adjustable Front Arms. The built-in rotating armrests will not increase the support width of the bed, but the rotating and swinging structure and locking mechanism will make removing the armrest operation easier, when doctors need space to reach the patient's chest and back, or to be closer to the patient.

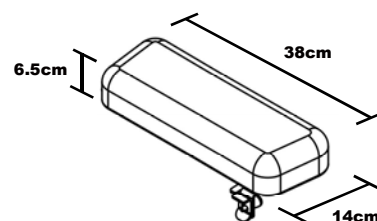

0 Without Anything

POSITION 8

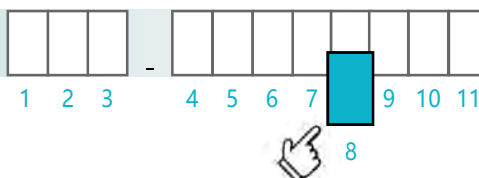

Stand-by Castor System




The liftable Stand-by Castor System adopts the medical castors and the connecting rod operating system. The four castors can be released and retracted synchronously by foot pedal control. The bed can be moved easily when wheels are released, and the stability of the bed will be increased when the castors are folded up. Castor system must be installed by the original factory. Users are not allowed to install the castor system without permission to avoid accidental injury.

1 With Stand-by Castor System

0 Without Stand-by Castor System

POSITION 9

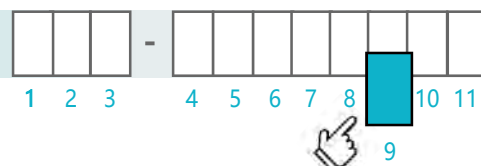

Electrical Safety Lock

The electrical safety lock adopts the traditional mechanical lock, which is safe and reliable. It is used to completely cut off the power supply of the motor and effectively avoid the operation and accidents of irrelevant personnel.




1 With Electrical Safety Lock

According to the European law, the electric lifting bed must be equipped with safety devices to avoid accidental injury caused by the operation of children or non-designated personnel. Each lockset is equipped with two keys.

0 Without

 The electrical safety lock must be installed by the original factory. It is forbidden for users to install it without permission. Please cut off the power supply and remove the key when the designated operator leaves.

POSITION 10

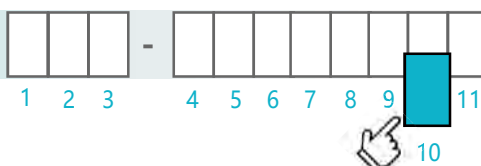

Paper Roll Holder




1 With Paper Roller Holder (Part code: PR02)

The Paper Roll Holder can be installed by users or factory. Maximum allowed width of paper roll is 60cm and maximum allowed diameter of paper roll 18cm, which you can purchase anywhere in the market. Hardware pack for installation is included.

 Paper Roll Holder is not applied on Treatment Avalon & Danvers.

0 Without

POSITION 11

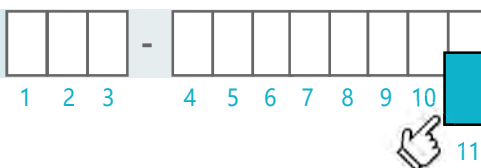

Sponge System



1 Soft Polyurethane Foam with PVC or PU Vinyl

Soft polyurethane sponge is usually used for massage, spa and as a home care bed. The thickness of the sponge pad is usually 2.5" or 3", which can use any leather Sintang provided. (Refer to Color Range page)

2 Hard Polyurethane Foam with PVC or PU Vinyl

Hard polyurethane sponge is usually used in clinics and professional treatment institutions. Sponge pads are usually 2" or 2.5" thick and can be made of any leather from Sintang. (refer to color range page)



3 Self-Skin Injection Foam

Our self skinning sponge pad is made of 200kg high-density polyether through high-pressure sealed injection molding, and the skin is integrally formed. It is featuring durable, waterproof, fireproof, oil resistance and anti-microbial. Self-skin sponge system is allowed for limited models and sizes. (Refer to Configuration Table)

 Self-Skin Injection Foam is only suitable for Treatment Cabell, Infinity, Scissors, Russell, Avalon, Danvers with width of 27" and 29".

Model Code	CTB	CTC	CTI	CTS
Image				
- Length	198cm			
- Width	22"/56cm, 25"/63cm, 27"/68cm, 29"/73cm			27"/68cm, 29"/73cm
- Product Weight	136lbs/62kg	169lbs/77kg	174lbs/79kg	176lbs/80kg
POSITION CODE 4 Control System				
A - Single Pedal	✓	✓	✗	✓
B - Double Pedal	✗	✗	✓	✗
C - Single Pedal + Side Panel Control	✗	✗	○	✗
D - Single Pedal + 2 Button Hand Control	○	○	✗	○
E - Single Pedal + 4 Button Hand Control	✗	✗	○	✗
F - Hand Control	○	○	○	○
G - Walk-Free Control	○	○	✗	○
H - Walk-Free Control + Side Panel Control	✗	✗	○	✗
I - Walk-Free Control + Hand Control	○	○	○	○
POSITION CODE 5 Breath Hole Design				
A - No Breath Hole	○	○	○	○
B - One Breath Hole + Filler	✓	✓	✓	✓
C - One Breath Hole + Neck Pillow	○	○	○	○
D - Two Breath Hole + Filler	✗	✗	✗	✗
E - Face Cushion + Filler	○	○	○	○
F - Two Breath Hole + Neck Pillow	✗	✗	✗	✗
POSITION CODE 6 Front Armrest				
1 - Height Adjustable Front Arms	○ 27", 29"	○ 27", 29"	○ 27", 29"	○ 27", 29"
2 - Rotatable & Height Adjustable Front Arms	○ 27", 29"	○ 27", 29"	✗	○ 27", 29"
POSITION CODE 7 Side Armrest & Guardrail				
1 - With Revolving Extended Side Arms	○ 27", 29"	○ 27", 29"	○ 27", 29"	○ 27", 29"
2 - With Revolving Guardrail	○ 27", 29"	○ 27", 29"	○ 27", 29"	○ 27", 29"
3 - With Rotatable Built-in Side Arms	○ 27", 29"	○ 27", 29"	○ 27", 29"	○ 27", 29"
POSITION CODE 8 Stand-by Castor System				
1 - With Stand-by Castor System	○	○	○	○
POSITION CODE 9 Electrical Safety Lock				
1 - With Electrical Safety Lock	✓	✓	✓	✓
POSITION CODE 10 Paper Roll Holder				
1 - With Paper Roll Holder	○	○	○	○
POSITION CODE 11 Sponge System				
1 - With Soft Polyurethane Foam	✗	✗	✗	✗
2 - With Hard Polyurethane Foam	✓	✓	✓	✓
3 - With Self-skin Injection Foam	✗	○ 27", 29"	○ 27", 29"	○ 27", 29"

☒ Basic Configuration
 ☐ Optional Configuration
 ☒ Unavailable

			-								
--	--	--	---	--	--	--	--	--	--	--	--

CTR	CTD	CTA	CBF	CBT
				
198cm			200cm	
22"/56cm, 25"/63cm, 27"/68cm, 29"/73cm			40"/100cm, 47"/120cm	
147lbs/67kg		167lbs/76kg	176lbs/80kg	185lbs/84kg
☒	☒	☒	☒	☒
☒	☒	☒	☒	☒
☒	☒	☐	☒	☒
☐	☐	☒	☐	☐
☒	☒	☐	☒	☒
☐	☐	☐	☐	☐
☐	☐	☒	☐	☐
☒	☒	☐	☒	☒
☐	☐	☐	☐	☐
☐	☐	☐	☐	☐
☒	☒	☒	☒	☒
☐	☐	☐	☐	☐
☒	☒	☒	☒	☒
☒	☒	☒	☒	☒
☒	☒	☒	☐	☐
☐	☐	☐	☒	☒
☐	☐	☐	☒	☒
☐	☐	☐	☐	☐
☐	☐	☐	☐	☐
☒	☒	☒	☒	☒
☐	☒	☐	☐	☐
☐	☐	☐	☒	☒
☒	☒	☒	☒	☒
☒	☒	☒	☒	☒

Camino Treatment Basic



CTB-AB000102



CTB-AE001102



CTB-GB130112



CTB-AB201112



CTB-FE221112



CTB-FB021102



CTB-GB231112



CTB-GE011112

Camino Treatment Cabell



CTC-AB000103



CTC-GE010113



CTC-AA020113



CTC-FB021103



CTC-AE200113



CTC-GB011113



CTC-GB111113



CTC-GE131103

Camino Treatment Infinity



CTI-BB000113



CTI-BE000103




CTI-CB110113



CTI-EE010113



CTI-FB121103

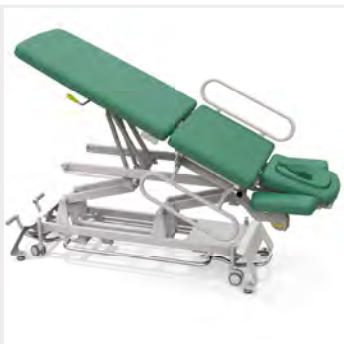

CTI-HE121103



CTI-IB011103



CTI-CE121113

Camino Treatment Scissors



CTS-AA000103



CTS-AB000103



CTS-FE010103



CTS-FB111103



CTS-GE230103



CTS-GB031103



CTS-GE021103



CTS-GB021103

Camino Treatment Russell



CTR-AB000102



CTR-FA010112



CTR-GB010102



CTR-GB121112



CTR-GB021112



CTR-GB211102



CTR-GB221112




CTR-FB110112

Camino Treatment Danvers



CTD-AB000103

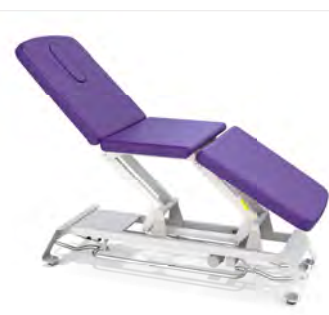

CTD-GB000103



CTD-GA010103



CTD-GB101103



CTD-GB121103



CTD-GA021103



CTD-GB211103



CTD-AB201103

Camino Treatment Avalon



CTA-BB000103

CTA-EB110103

CTA-CA000103

CTA-FB021103



CTA-HA021103

CTA-IB211103

CTA-HB011103

CTA-BB020103

Camino Bobath Flat & Bobath Tilt



CBF-AA000112

CBF-FB000112

CBF-GD021112

CBF-GF001112



CBT-GA001102

CBT-GB021112

CBT-GF021112

CBT-FD01112

CAMINO MASSAGE TABLES



Breast Comfort System




Paper Roll Holder



Stand-by Castor System



Independent Bolster System

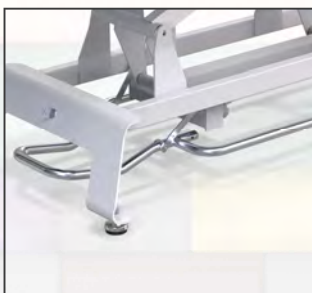

Walk-Free Height Control



Electrical Safety Lock

Camino Massage Flat

CMF -


--	--	--	--

The Camino massage table meets the concerns and needs of serious massage therapists and high-end spa club owners, who want to provide their customers the best in comfort, support and upscale atmosphere. It uses an exclusive segmented combination mattress that provides the comfort of a high-end couch and provides better breathability for the contact surface. Split mattresses can be arbitrarily combined according to user needs. Exclusive Chest Recesses and Neck Contour features provide your male and female clients with unparalleled comfort during face-down and face-up massages.



Specification

- Length: 72.5"/184cm
- Width: 28"/69cm, 30"/74cm
- Height: 18"/46cm~37"/95cm (Dewert)
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 167lbs/76kg

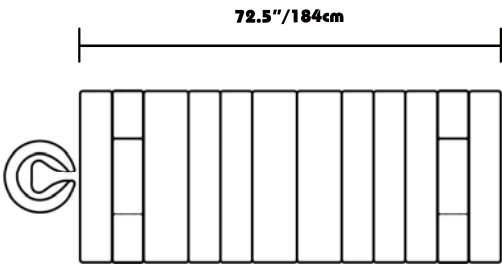

Design Varieties:



CMF-A010



CMF-F010



CMF-G110



CMF-G010



CMF-G011



CMF-G111



CMF-F011



CMF-A011

Camino Massage Tilt




CMT is a multi-purpose bed, suitable for massage, beauty and SPA use. CMT not only meets the concerns and needs of serious massage therapists and high-end spa club owners, but also adds a practical pneumatic backrest feature. Its unique 4" / 10cm thick sectional combination mattress can provide the comfort of a high-end sofa alike and provide better breathability for contact surfaces. The sectional mattress is composed of 16 separate removable and replaceable round pillows, which can be arbitrarily combined according to user needs. Exclusive Chest Recesses feature provides both your male and female clients with unparalleled comfort during face-down and face-up massages.



Specification

- Length: 72.5"/184cm
- Width: 28"/69cm, 30"/74cm
- Height: 18"/46cm~37"/95cm (Dewert)
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 167lbs/76kg

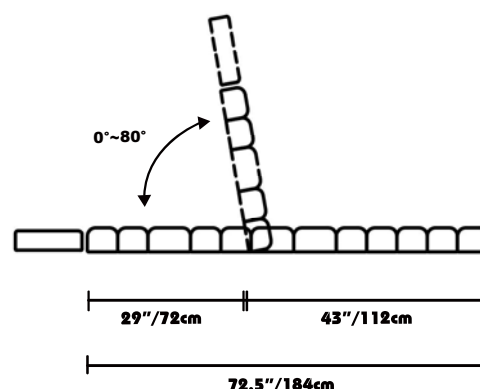

Design Varieties:



CMT-G010



CMT-A010



CMT-F010



CMT-G110



CMT-F011



CMT-A011



CMT- G011



CMT-G111

Camino Massage Deluxe


C	M	D	-				
1	2	3	4	5	6	7	

CMD is a high-end electric SPA massage bed. The 4-section adjustable mattress can provide users with a zero-gravity comfortable sitting position. The table mattress adjustment ability can meet the treatment needs of different supine or prone directions. Its exclusive 4" / 10cm thick sectional combination mattress can provide the comfort of a high-end sofa alike and provide better breathability for contact surfaces. The sectional mattress is composed of 16 separate removable and replaceable round pillows, which can be arbitrarily combined according to user needs. With breast recesses and padded pads at both ends of the bed, it offers unparalleled comfort for both male and female customers.



Specification

- Length: 72.5"/184cm
- Width: 28"/69cm, 30"/74cm
- Height: 18"/46cm~37"/95cm (Dewert)
- Foam: 2.5"/6.4cm
- Lifting Capacity: 440lbs/ 200kg
- Working Weight: 600lbs/272kg
- Product Weight: 167lbs/76kg

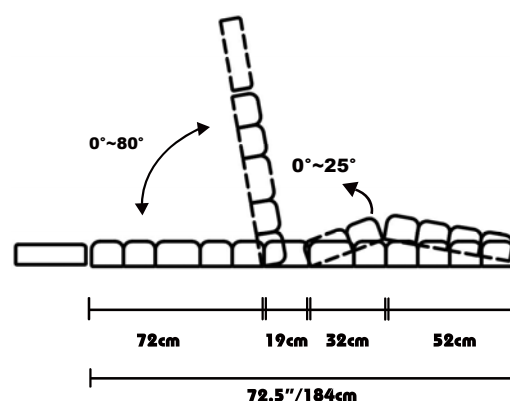

Design Varieties:



CMD-B010



CMD-F010



CMD-I110



CMD-H110



CMD-F011



CMD-B011



CMF-H110



CMF-I010

Meaning of Product Code

1	2	3	4	5	6	7	



Model Code	CMF	CMT	CMD
Image			
POSITION CODE 4 Control System			
A - Single Pedal	✓	✓	✗
B - Double Pedal	✗	✗	✓
C - Single Pedal + Side Panel Control	✗	✗	○
D - Single Pedal + 2 Button Hand Control	○	○	✗
E - Single Pedal + 4 Button Hand Control	✗	✗	○
F - Hand Control	○	○	○
G - Walk-Free Control	○	○	✗
H - Walk-Free Control + Side Panel Control	✗	✗	○
I - Walk-Free Control + Hand Control	○	○	○
POSITION CODE 5 Stand-by Castor System			
1 - With Stand-by Castor System	○	○	○
POSITION CODE 6 Electrical Safety Lock			
1 - With Electrical Safety Lock	✓	✓	✓
POSITION CODE 7 Paper Roll Holder			
1 - With Paper Roll Holder	○	○	○

Camino Mattress Bolster Cube

Camino's exclusive segmented combination mattress can be combined according to users' requirement. Each bolster is an independent individual. The bolster uses hard sponge as the pillow core to improve body support. The soft mercerized filler wrapped in the outer layer improves the comfort. Bolsters can be sold separately, 10 pieces / package.



#1 Camino Stripe Bolster Cube

19cm * 10cm * 76cm (Model: CMB#1-30)

Suitable for 30" Camino massage table.
Available in multiple colors. MOQ: 10 pieces

19cm * 10cm * 70cm (Model: CMB#1-28)

Suitable for 28" Camino massage table.
Available in multiple colors. MOQ: 10 pieces



#2 Camino Stripe Bolster Cube

13cm * 10cm * 76cm (Model: CMB#1-30)

Suitable for 30" Camino massage table.
Available in multiple colors. MOQ: 10 pieces

13cm * 10cm * 70cm (Model: CMB#1-28)

Suitable for 28" Camino massage table.
Available in multiple colors. MOQ: 10 pieces



#3 Camino Small Bolster Cube

13cm * 10cm * 22cm (Model: CMB#3-30)

Suitable for 29" Camino massage Table.
Available in multiple colors. MOQ: 10 pieces

13cm * 10cm * 19cm (Model: CMB#3-28)

Suitable for 27" Camino massage Table.
Available in multiple colors. MOQ: 10 pieces



#4 Breast Recesses Filler

13cm * 10cm * 32cm (Model: CMB#4)


Suitable for any size of Camino massage Table. Available in multiple colors. MOQ: 10 pieces



#5 Breast Recesses Padded Pad




13cm * 4cm * 32cm (Model: CMB#5)

When massaging female with smaller breasts, fill 1-2 sponge pads in the breast recesses.
Available in multiple colors. MOQ: 10 pieces

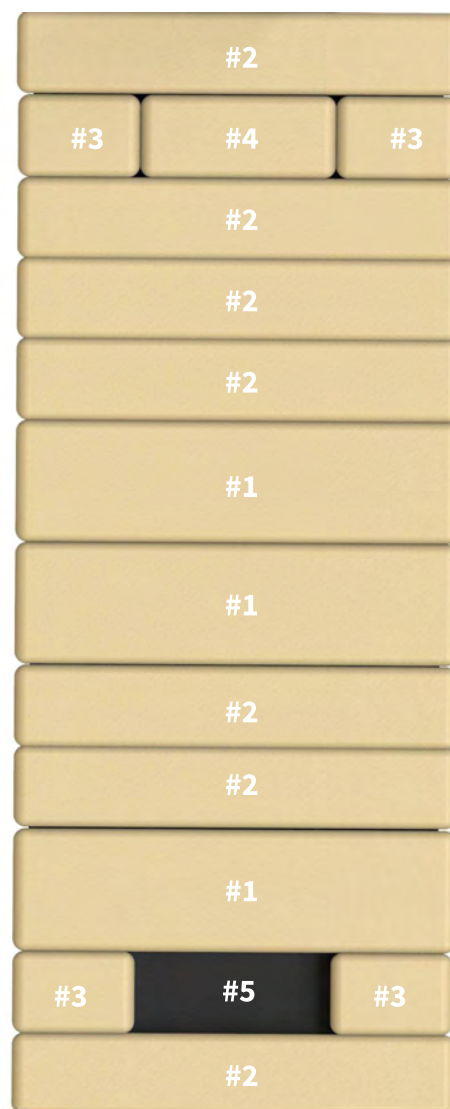

● Standard PU ● Eco PU ● Italino PU ● PVC

Standard PU leather supports 5-year hydrolysis resistance; Eco Pu supports 2-year hydrolysis resistance; Italiano high gloss PU supports 5-year hydrolysis resistance. Although PVC material is not environment friendly, it can be used for longer period of time.

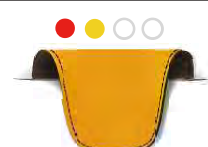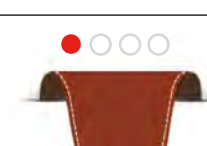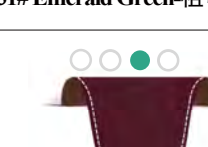 01# Black-黑色	 02# Burgundy-葡萄红	 03# Purple-紫色	 04# Agate Blue-灰色	 05# Teal-鸭绿
 06# Beige-米色	 07# Mc Beige-米黄	 08# Chocolate-巧克力	 09# Navy Blue-海军蓝	 10# Lilac-粉红色
 11# Blue-蓝色	 12# Yellow-黄色	 14# White-白色	 15# Macaw Green-棕榈绿	 16# Jade-玉色
 17# Sky Blue-天蓝	 18# Seamist-浅蓝灰	 19# Orange-橙色	 20# Cyan-青色	 21# Otter-水獭色
 22# Charcoal-炭灰色	 23# Coffee-亮咖	 24# Cinnamon-肉桂色	 25# Cobalt Blue-艳蓝色	 26# Lily Green-淡绿色
 27# Surf Blue-水蓝	 28# Dove Grey-鸽子灰	 29# Cream-亮米黄	 30# King's Purple-亮紫色	 31# Emerald Green-祖母绿
 32# Maina Blue-玛丽娜蓝	 33# Formula One-赛车红	 34# Ochve-黄褐色	 35# Meadow Green-草绿	 36# Dark Purple-深紫
 37# Tobacco Brown-烟棕色	 38# Maroon-褐红色	 39# Cabernet-胭脂色		

No.136 Linsheng Road, Jinshan District, Shanghai



Factory:

Tel: **+86-21-37210909**

Add: **No.136 Linsheng Road, Jinshan District,**
Shanghai, China



Show Room U.S.A:

Tel: **1-888-710-7206**

Add: **8360 Rovana Circle, Sacramento,**
California, 95828, U.S.A



Show Room UK:

Tel: **+44 (0)772 1762 663**

Add: **Unit 8 Fareham Industrial Park,**
Standard Way, Fareham PO16 8XB
United Kingdom



Show Room Germany:

Tel: **+49 (0) 151 2435 0568**

Add: **Hittfelder Kirchweg 21, 21220**
Seevetal, Germany



HEMET
— 禾玛 —

www.mingtenghealth.com