


# **“Patented” Turning Mattress**

## for Anti-Bed Sore & Pressure Ulcer Prevention

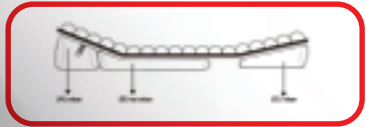


CE FDA ISO 13485

## Rotation mattress for disabled patients



**HF808**

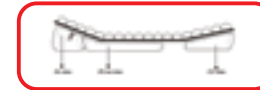

### Features and Benefits

- Tubular made of TPU film/Nylon with TPU coating
- AB alternative for ulcer prevention
- Welded, seam sealed, zipper around waterproof and vapour permeable stretch PU cover
- Easy turn CPR for rapid deflation
- Silent running pump at optimum support pressures
- Optional 12/18/24 mins cycle time
- Visual and audible alarms for patient safety
- With remote control
- External, Easy-to-replace pump filters with static function and auto-return
- Automatically adjusts to the patients weight with head up, foot up, left rotate, and right rotate functions

Item:HF808

- Pump Size:28(L)\*10(W)\*16.5(H)cm
- Minimum Weight Limit:No Minimum
- Maximum Weight Limit:135-145kg
- Mattress Size:200\*90\*11.5-40cm
- Material:TPU film / Nylon TPU for air cells,Stretch PU
- Cover,Bottom Tarpaulin/Nylon PU

